



Attention!

Please read this document
before unpacking your

HunterLab[®] **ColorFlex[®] L2**

A60-1021-432, ver 1.0



1. Unpack Your Box

Place the ColorFlex L2 on the bench. Remove wrappings and cable ties. Retain the packaging for return of the instrument to HunterLab if needed.

2. Review Your Order

Please review your order. Compare your purchase order with the shipping document from HunterLab. If there are any discrepancies, please contact the HunterLab Customer Experience Team at <https://www.hunterlab.com/en/contact/>.


3. Select a Space for the ColorFlex L2

Laboratory grade environment is required to ensure precise and accurate measurements. The table below shows instrument size and weight, electrical, and temperature/humidity requirements.

| | |
|--------------|---------------------------|
| Height | 16 cm; 6.25 in |
| Width | 21 cm; 8.75 in |
| Depth | 24 cm; 9.5 in |
| Weight | 2.75 kg; 6 lb |
| Power Outlet | 110-240AC |
| Electrical | 47-63Hz Single Phase |
| Temperature | 10 - 30°C (50 - 86°F) |
| Humidity | 10% to 90% non-condensing |


4. Site Conditions

- Allow 20 cm (8 in) of unobstructed clearance on each side of the instrument.
- The Instrument should be placed to **avoid direct sunlight** and away from windows.
- The atmosphere around the ColorFlex L2 should be **free from alkalis and acid vapors, and organic solvents.**
- Do not operate the ColorFlex L2 in the presence of **strong electrical fields.**
- **Avoid strong vibrations or shocks.**
- Ensure that there are no **heat-generating sources** near the ColorFlex L2.

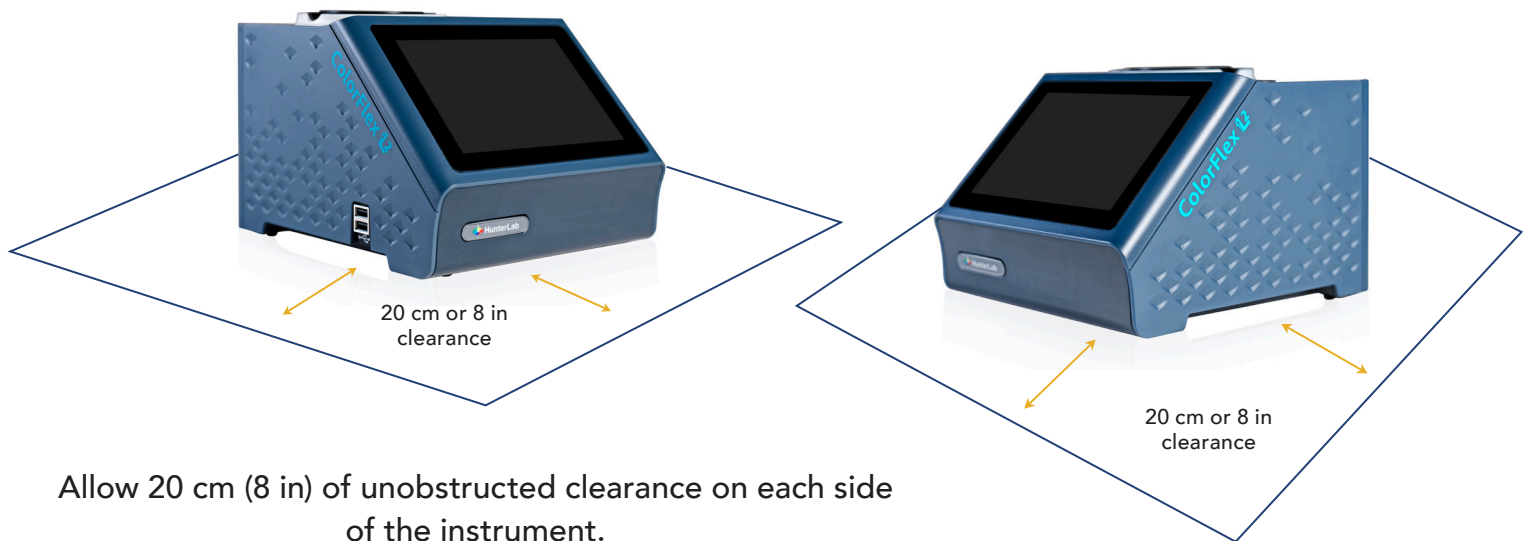

Allow 20 cm (8 in) of unobstructed clearance on each side of the instrument.

Note: Failure to comply with these conditions and protocols set forth in this document may adversely affect the instrument performance.

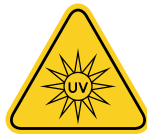

5. Safety Requirements

- UV Radiation – Wear UV protecting eye wear when viewing or observing the LED's inside the instrument. The frequency of this flashing light is in the range of sensitivity for those prone to epileptic seizures.

6. Instructions

HunterLab has included a printed Quick Start Guide with the shipment. This guide is intended to get you started with measurements. Also included is the HunterLab Instruction Manual for ColorFlex L2 on a flash drive. This manual provides a detailed view of operation, troubleshooting and optimization. By referring to the instruction manual, you will gain understanding of practical applications that will help you achieve excellence in using your HunterLab ColorFlex L2 spectrophotometer.

7. Context Sensitive Help

The context sensitive help for ColorFlex L2 complements the information presented here and supplies a more detailed view of operation by accessing information screen by screen. By referring to the on-board help, you will gain understanding of practical applications that will help you achieve excellence in using your HunterLab ColorFlex L2 spectrophotometer.

8. Service

- Only HunterLab authorized, factory-trained representatives are permitted to service HunterLab's products.
- Do not remove any of the covers. Service is to be performed by a HunterLab service representative.

Contact the HunterLab Customer Experience Team at support@hunterlab.com if service is needed.

For more information on the ColorFlex L2, scan the QR code.



SCAN ME



11491 Sunset Hills Road
Reston, VA 20190-5280
www.hunterlab.com