

D25[®] NC
FEATURES



Sensor



Sealed aluminum case and optics to protect sensor in dusty environments.



LED Light source
is bright and
reliable
(5 yrs warranty.)



HunterLab



Laser Distance
Sensor - Product
Height or Distance
from sensor data to
ensure correct
sample presentation.

 HunterLab



Enter Password



DONE +
-



HunterLab



250 stored
Standards with
password
protection

Saved Readings

DELETE ALL READINGS ▶

PRINT ALL READINGS ▶

VIEW SAVED READINGS ▶

FILTER BY SETUP ▶

*Datalog holds
previous 1000 color
data measurements
for onscreen review
or digital download*



HunterLab



Standardization


| | |
|------|-----------------|
| 0.05 | PASS |
| 0.03 | |
| 0.06 | (≤ 0.30) |

Final Result

BACK  MAIN MENU 

Keypad and display
for standardization,
color measurement,
and system
diagnostics.

 **HunterLab**

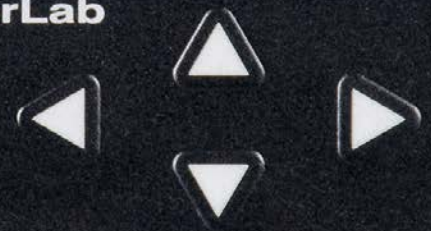

Data View

| | |
|----------------|----------|
| VIEW 3 | ENABLED |
| DISPLAY | ABSOLUTE |
| ILL/OBS | C/2 |
| COLOR SCALE | L a b |
| COLOR INDEX | dE |
| CMC/SHADE SORT | ▶ |



HunterLab



Powerful color calculator for a variety of color scales, indices, and illuminants.

D25 NC Base



5 seconds

*Complete
rotation in
5 seconds.*

D25 **NC**


*Automatic
Start/Stop
rotation for each
scan, with motion
disable switch.*

MOTION

OFF



AUTO


Sealed transport
mechanism.



*Hand disassembly
for easy cleanup
of spills.*



| GREEN | WHITE |
|----------------------|---------------------|
| STANDARD | STANDARD |
| Aug 11, 2016 | Aug 21, 2016 |
| X = 14.70 | X = 79.61 |
| Y = 20.63 | Y = 84.56 |
| Z = 16.63 | Z = 86.83 |
| Diagnostic: Standard | Redundant: Standard |



Direct drive
motor for
samples weighing
up to 20 lbs.



USB port for
connection to
a printer, barcode
scanner. Compatible
with EasyMatch QC
4.8 and above on
Windows based PC.



MAIN POWER

*A/C power
switch with
detachable
International
power cord
(IEEC.)*



100-230
50/60 HZ



3AG 2A
X2

*Hardened
powder-coated
aluminum chasis
for stability and
durability.*

D25 NC


MOTION
OFF  AUTO



Adjustable Support Arms




Height Adjustment
to accommodate
samples from 2.5
to 5.5 inches from
the sensor window.



*Width Adjustment
for a range of
sample dish sizes
from 2 to 12 inches.*

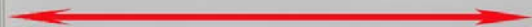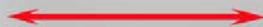

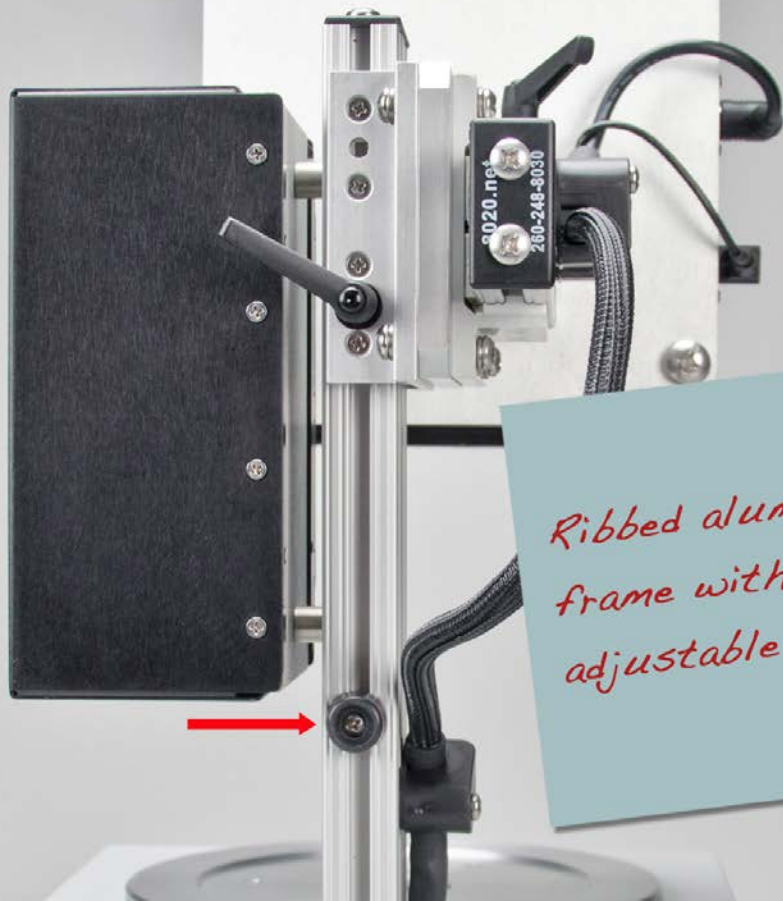

*Ribbed aluminum
frame with
adjustable stops*

Sample Dishes

D25 **NC**


*Includes 6-inch
sample dish*





D25 NC

*Includes 12-inch
sample dish*

Sample dishes can
be nested to use
maximum scan area
with minimal sample.



Utility Software

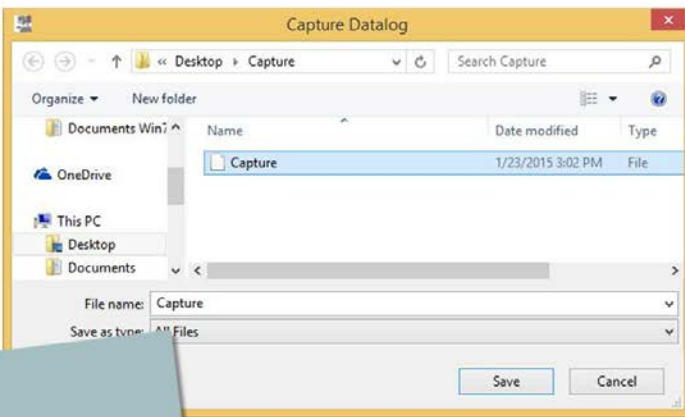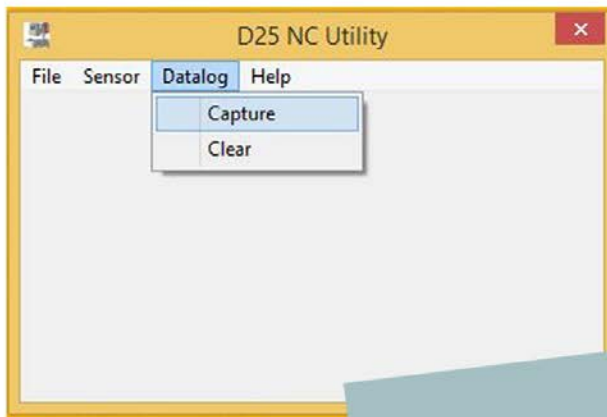

Retrieve sample measurement datalog.

D25 NC Utility

File Sensor Datalog Help

Archive Product Setups

Restore Product Setups

Edit Product Setup

Identification

Set Date and Time...

Archiving Product Setups...



Cancel

69 / 250

Archive Product Setups

Desktop > Capture

Organize New folder

| Name | Date modified | Type |
|-----------------------------|---------------|------|
| No items match your search. | | |

File name: Capture

Save as type: All Files

Save

Cancel

*Backup storage
and retrieval
of Product
Standards.*

D25 NC Utility

File Sensor Datalog Help

Archive Product Setups

Restore Product Setups

Edit Product Setup

Identification

Set Date and Time...

Edit Product - D25 NC

Name: SETUP 1

Standard: HITCH

Autosearch:

Views: 1 Edit

Standard Values: Edit

Tolerances: Edit

Run Configuration: Edit

1

Load

Save

Cancel

Product Standards editing capability.

For more information go to

www.hunterlab.com

or contact your local HunterLab representative.